

MediCap USB200



Ref.-No. 190.0020


MediCapture USB200, Medizinisches digitales Bild- und Videospeicher-Gerät mit USB-Funktion
 Einsetzbar mit nahezu allen gängigen Videobildgeräten (in Endoskopie, Dental, Veterinär, Ultraschall, ...)

Von analogen medizinischen Bildgeräten werden Bilder und Videos per Knopfdruck oder über Fußschalter digital erfasst und gespeichert. Ausgewählte Bilder und Videos können über den integrierten USB-Port an einen USB-Stick gesandt und dort gespeichert werden. Dieser Stick ist zum Anschluss an alle handelsüblichen USB-kompatiblen Geräte wie z. Bsp. PC, Laptop, usw. geeignet. Über die im Stick integrierte einfache Software können die Bilder und Videos dann bequem am jeweiligen Gerät (PC, Laptop, ...) verwaltet werden.

MediCapture USB200, Digital Medical Image and Video Capture Device with USB-Function
 For use with nearly all usual video imaging devices (in Endoscopy, Dental, Veterinary, Ultrasound, ...)

From analog medical-image-devices digital images and videos can be easily captured and saved by button or over footswitch. Chosen images and videos can be transferred through the integrated USB-Port to the USB-Flash-Drive-Stick and additionally can be saved there also. This stick is suitable for the connection to all usual USB-compatible devices like PC, Laptop, ... ! Comfortable storage of the image and video files at your PC, Laptop, ... is possible through the stick-integrated easy software.

Rückseite /
 Back panel :

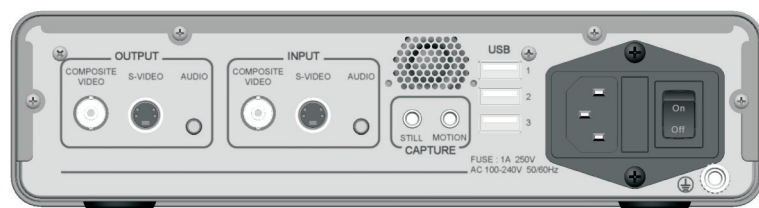

Technische Daten / Technical Data:

2,5" Monitor integriert /integrated (for navigation only)	Abm. (BxHxT) / Dim. (WxHxD):	240 mm x 63 mm x 200 mm ; 1,5 kg
TV-System: PAL/NTSC - umschaltbar - switchable	Temperaturen /	Lagerung/Storage: - 40° / +85° C
Dateiformate / File formats: JPEG, TIFF, PNG / MPEG2, MPEG1	Environmental:	Betrieb/Operating: -20° / + 40° C
Bildgröße / Image size: 1024 x 768, 800 x 600, 640 x 480 Pixel	Lieferumfang/Contents:	MediCapture, USB-Stick (1 GB)
Kamera Input / Output: S-Video, BNC	Mit integrierter/ with integrated Software	UL, FCC, CE, FDA Class 1,
Netzspannung / Power supply: 100 VAC - 240 VAC, 50-60 Hz	Zertifikate / Certificates:	EN60601-1, EN60601-1-2
	optional:	Fußschalter / Footswitch 190.0000.FS

Jetzt auf Anfrage auch erhältlich:

Now also available on request :

MediCap USB300 HD (High Definition), Ref.-No. 190.0020.HD